ORIGIAL

ORICAREFIT-W

ORIGINAL SCALE CARE&FIT SMART WHITE

USER MANUAL

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EN

EN

ORIGINAL SCALE CARE&FIT SMART WHITE OPERATING AND SAFETY INSTRUCTIONS MODEL NO. ORICAREFIT-W

FOR INDOOR HOUSEHOLD USE ONLY

SAVE THESE INSTRUCTIONS

This appliance is not intended to be used in household and similar applications such as:

- staff kitchen areas in shops, offices and other working environments;
- farm houses;
- by clients in hotels, motels and other residential type environments;
- bed and breakfast type environments.

CONTENT

Introduction	4
Product Overview	4
Technical Specifications	5
Installation	5
Startup	5
Using the Scale	6
Safety Warnings	6
Scale Maintenance	8
Warranty	8

**PLEASE READ THIS MANUAL CAREFULLY
BEFORE YOU USE THE APPLIANCE AND SAVE
IT FOR FUTURE REFERENCE.**

Introduction


The Original Scale Care&Fit Smart is a smart digital scale designed for comprehensive monitoring of your weight and body composition at home. Thanks to its high-precision sensors and Bluetooth connectivity, it allows you to measure not only body weight but also parameters such as fat percentage, muscle mass, and other biometric values, synchronizing them with your phone through the Feelfit app for detailed health tracking. This premium design scale combines a durable tempered glass platform with stainless steel electrodes to deliver reliable results and a modern look. It is an easy-to-use tool, designed for the whole family, that will help you monitor your well-being comfortably and safely every day.

Product Overview

Figure 1: Top view of the scale, showing the tempered glass platform, measuring electrodes, and illuminated LED display.

The main parts of the Care&Fit Smart scale are as follows:

- 1.- Measuring surface: Tempered glass weighing platform where the user stands.
- 2.- Measurement contact area: Metal electrodes that must be touched with bare feet for body composition analysis.
- 3.- LED Display: Shows weight and status symbols (e.g., Bluetooth indicator, overload, low battery, etc.).
- 4.- Reset button: Located at the bottom, used to reset or calibrate the scale if necessary.
- 5.- Battery compartment: Space at the base of the scale to install the batteries that power the device.
- 6,7.- Non-slip feet: Support legs with rubber material providing stability and floor protection, housing the weight sensors inside.
- 6,7.- Weight sensors: Internal high-precision components (strain gauge type) located in the feet, responsible for measuring load.
- 8.- Bluetooth indicator: Icon on the screen that lights up when the scale is properly connected via Bluetooth to the mobile device.

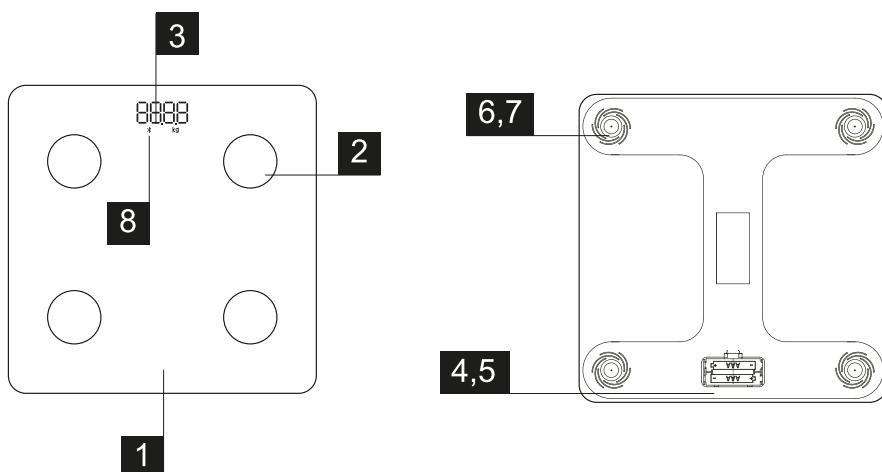APP FEELFIT

EN

Technical Specifications

Battery type	2 × 1.5 V AAA batteries (placement in bottom compartment)
Display	High-visibility LED display with red digit illumination
Sensors	4 high-precision weight sensors (strain gauge technology)
Weight capacity	Maximum weight 180 kg (396 lb)
Measurement unit	Kilograms (kg) by default; switchable to pounds (lb) via the Feelfit app
Wireless connectivity	Low-power Bluetooth® 4.0 (BLE) for smartphone synchronization
Compatible app	Feelfit (available for iOS and Android, free download)
Materials	Durable tempered glass platform; stainless steel electrodes; ABS plastic base and feet
Auto on/o	Yes – the scale turns on automatically when weight is detected and turns off after a few seconds without use
Other functions	BMI calculation and body composition analysis in the app; storage of multiple user profiles
Note	Weight accuracy is ±0.1 kg. This scale is designed for home use only. Body composition data provided are for reference and do not replace professional medical advice.

Installation

Follow these steps for correct installation and preparation of the scale before first use:

- Battery placement: Open the battery compartment cover (5) at the base of the scale and insert 2 new AAA batteries, respecting the indicated polarity. Make sure they are properly fitted.
- Proper location: Place the scale on a flat, solid, and level surface. Avoid placing it on thick carpets, uneven or slippery surfaces, as these could affect measurement accuracy.
- Initial check: Before using, and with the scale in a horizontal position, step briefly on the surface (1) and step off to activate and calibrate the device. The display should turn on and then show "0.0" (indicating it is ready to measure). If the scale does not turn on, check battery installation or gently press the reset button (4) with a thin object.

Startup

Once installed, proceed to set up and link the scale with your smartphone:

App download: On your mobile phone, download and install the Feelfit application. It is available for free in the iOS App Store and Google Play for Android. You can also scan the QR code included in the product documentation to directly access the download.

Bluetooth activation: Turn on Bluetooth connectivity on your smartphone. Open the Feelfit app and follow the initial instructions to register a user or log in if necessary.

Pairing with the scale: With the app open near the device, step briefly on the scale to turn it on. The scale will enter synchronization mode and display the Bluetooth symbol (8) on the screen. The scale should appear in the app; select the ORICAREFIT-W model or the equivalent option to pair. The connection will be established in a few seconds. Note: The scale may require you to step on it during initial pairing to activate.

EN

Confirmation: Once paired, the Care&Fit Smart will remain connected to the phone during use. The Bluetooth connection icon will light up on the display (3) whenever the scale transmits data to the app. Pairing is only required the first time; in later uses, the connection will be automatic. If you ever need to link it to another device or solve connection issues, hold the reset button (4) for a few seconds to reset Bluetooth settings and repeat the pairing process.

Using the Scale

To obtain accurate measurements of weight and body composition, follow these usage guidelines:

- User placement: Step on the scale barefoot, placing both feet on the measuring surface (1) and ensuring each foot makes good contact with the metal electrodes (2). Stand in the center of the platform, with feet evenly apart. It is important that the skin of your feet is clean and slightly moist (not excessively dry) to improve conductivity for bioimpedance measurement.
- Posture and waiting: Remain upright and completely still while the scale performs the measurement. The LED screen (3) will first display body weight and then an animation or progress indicator for a few seconds while calculating other parameters (fat, muscle, etc.). Avoid moving or speaking until the measurement is complete.
- Result reading: Once the measurement is complete, the screen will display the final weight. Other detailed body composition data will automatically sync with the Feelfit app on your smartphone via Bluetooth (the app must be open or running in the background with an active session). In the app, you can consult values such as body fat percentage, BMI, muscle mass, water percentage, among others, with progress charts and recommendations.
- Completion: After stepping off the scale, it will turn off automatically after a few seconds of inactivity to preserve battery life. Manual shutoff is not necessary. Simply store the scale in a safe, dry place until the next use.
- Usage tips: For consistent results, try to weigh yourself at the same time of day (e.g., in the morning) under similar conditions. Always place the scale on a hard floor (avoid carpet or soft surfaces) to ensure accuracy. If the scale has been moved or stored vertically, it is recommended to perform a quick calibration by stepping on and off the platform before the next measurement. Keep the electrode surfaces clean with a dry cloth to ensure good conductivity (do not use chemicals or immerse the device in water).

Safety Warnings

Please read carefully the following safety warnings and use them as a reference throughout the product's life. Failure to follow these precautions could result in injury to the user or damage to the scale.

Danger! People with pacemakers or other implanted electronic medical devices must not use the body composition analysis functions of this scale. The bioimpedance technology sends a very low electrical current through the body that could interfere with such devices. Consult your doctor before using it solely for weight measurement, or refrain from using this product.

Danger! Keep the scale's batteries out of reach of small children. Swallowing batteries can cause severe or fatal injuries. If you suspect a child has swallowed a battery, seek medical

EN

attention immediately. Also, ensure proper battery installation and never throw them into fire (risk of explosion).

Warning: This scale is designed for indoor household use only. Do not use outdoors or in very humid environments (e.g., bathrooms with excessive steam) to avoid electronic failures or short circuits. Do not immerse the scale in water or allow liquid inside the device. For cleaning, follow the maintenance instructions.

Warning: Do not place heavy objects on the scale when not in use, and never exceed the maximum capacity of 180 kg, as this could damage the weight sensors or break the glass surface. If the display shows an overload error message (e.g., "Err"), remove the weight immediately.

Warning: Use the scale barefoot for bioimpedance measurements. Do not weigh yourself wearing shoes or socks, as the electrodes will not obtain a reading and body composition data will be incorrect. Also, do not step on the scale with wet feet or when the surface is wet; it may be slippery and pose a fall risk. Dry both the platform and your feet before use.

Caution: The glass platform is designed to support the intended weight but can break if subjected to strong impact (e.g., if the scale is dropped). Handle the device with care. If you notice any cracks or damage in the glass, stop using the scale immediately to avoid cuts or injuries, and contact technical service for replacement.

Caution: Place your feet centered and balanced on the scale. Avoid standing on the edge or in an unstable position. Always use the scale on level surfaces; if placed on an uneven surface, there is a risk of tipping or inaccurate reading.

Caution: Do not jump or bounce on the scale, nor use it as a platform or step. Misuse may damage the sensors or cause accidents.

Caution: This device is not intended for use by very young children or people with disabilities without proper adult supervision. Teach older children how to use the scale correctly and safely, and prevent them from playing with it.

Important: Do not attempt to open or disassemble the scale. It contains no user-serviceable parts (except for the battery cover). Any repair must only be carried out by authorized technical personnel. Altering or modifying the product will void the warranty.

Important: If the scale will not be used for an extended period, remove the batteries to avoid damage from possible leakage. Store the device in a dry, safe place, away from heat sources or prolonged direct sunlight.

Important: This scale is a personal wellness assessment device. It is not a medical device. The body composition data provided are estimates based on population averages. Consider these results as guidance and consult a healthcare professional to interpret your particular situation.

EN

Scale Maintenance

To ensure the correct functioning of your Care&Fit Smart scale over time, please take into account the following maintenance recommendations:

- Battery replacement: When the low battery indicator appears on the screen (for example, an empty flashing battery symbol) or if you notice the screen illumination is dim, replace all the batteries with new ones of the same type (AAA). Do not mix used and new batteries, nor batteries of different brands, to avoid possible leaks or damage. Always insert the batteries respecting the indicated polarity.
- Storage and conservation: After each use, store the scale in a dry, clean, and safe place. Avoid storing it in a vertical position for long periods (ideally, keep it flat) and do not place heavy objects on it while stored. Protect it from extreme humidity and from very high or very low temperatures that could affect its electronic components.
- Automatic calibration: The scale is designed to calibrate automatically. If you have moved or transported it, it is recommended to reset it before the next use. To do this, with the scale on the floor, press the reset button (4) on the underside using a pointed object (for example, a paperclip) or simply step on and off the scale quickly to reset it to zero. When doing so, the display should show "0.0" indicating that it is ready and calibrated.
- General cleaning: (Note: No special frequent cleaning is required beyond the basics.) Keep the surface free of dust and dirt by gently wiping with a dry or slightly damp cloth. Ensure that no water enters the electrode holes or the battery compartment. Do not use solvents or abrasives that could damage the glass coating or plastic parts.
- Periodic inspections: Occasionally check that the non-slip feet (6) are clean and in good condition. If any become loose or lose grip, it could affect stability and measurement; adjust or replace them through technical service. Also, verify that the battery compartment cover closes correctly.

By following these guidelines, your scale will remain in optimal condition, ensuring precise measurements and extending its useful life. If you have any doubts or unresolved problems despite these precautions, please contact our technical support service.

Warranty

Terms and Conditions of the Warranty: This product comes with a limited 3-year warranty starting from the date of the consumer's original purchase. The warranty covers manufacturing defects in materials and workmanship under normal household use conditions. MLTEC MOTION INNOVATION S.L.U., as the authorized importer and distributor, certifies that the Original Scale Care&Fit Smart (ORICAREFIT-W) scale is free of such defects at the time of purchase. If, during the warranty period, the product exhibits failures due to defective materials or incorrect manufacturing processes, the consumer will be entitled to free repair or replacement of the product, under the conditions detailed below. The battery included with the product has a specific 6-month warranty, as it is a consumable component. This warranty is additional and does not affect the consumer's rights established by national legislation.

Warranty Coverage:

EN

The warranty will only be valid upon presentation of the stamped Warranty Certificate (or the original proof of purchase) together with the invoice or purchase receipt indicating the model purchased and the purchase date. The certificate or proof of purchase must clearly include: (a) the buyer's name, (b) the name, address, and stamp of the selling establishment, (c) the product model, and (d) the purchase date. The absence or alteration of any of this information may void the warranty service.

The following cases are expressly excluded from this warranty:

- Products whose warranty period (3 years) has expired.
- Damage or malfunctions caused by improper use or contrary to the instructions in the user manual (for example, overload above the maximum weight, exposure to liquids, falls, impacts, etc.).
- Unauthorized interventions, repairs, or modifications: any handling carried out by technicians not authorized by MLTEC MOTION INNOVATION, the installation of non-original parts or software, or alteration of the product's original specifications will void the warranty.
- Natural wear and tear from normal use of the product in consumable or perishable components. For example, elements such as depleted batteries (after 6 months), worn rubber or plastic parts due to regular use, scratches, or aesthetic discoloration over time are not covered.

Accidental damage or damage caused by external events: floods, fires, power surges, lightning, natural disasters, or other force majeure causes not attributable to product manufacturing defects.

If the product shows a defect covered by the warranty, the procedure to follow is as follows: the consumer must notify MLTEC Technical Support Service (see Technical Support section) within a reasonable time (preferably within 15 days from detection of the issue). Instructions will be provided on how to send the defective product, properly packaged, to the technical service. After evaluation, MLTEC may choose to repair the product (e.g., by replacing damaged components with new or refurbished parts) or replace it with a product of equal or equivalent performance. The chosen option will be free of charge for the consumer, including labor and materials.

Any repair or replacement carried out under warranty will only cover the remaining time of the original warranty period. That is, a new 3-year period does not begin; the period continues from the initial purchase date. All replaced parts and products will become the property of MLTEC MOTION INNOVATION.

Territorial Scope: This warranty applies only within the territory of the Kingdom of Spain. Warranty services will be provided within this territory. For products purchased outside Spain or if the user is located in another country, they must contact the corresponding local distributor to know the warranty coverage applicable in their region.

For any questions regarding the terms of this warranty, you may contact our customer service. This warranty acknowledges the consumer's legal rights under current regulations, without restricting them.

EN

Technical Support Service: If you need assistance with the product's operation, to make the warranty effective, or to purchase original spare parts, you can contact our customer service and technical support through the following means:

Email: soporte@tec1.zendesk.com

Website: www.unotec.es
(Support/Contact section)

ES

ORIGINAL SCALE CARE&FIT SMART WHITE OPERATING AND SAFETY INSTRUCTIONS MODEL NO. ORICAREFIT-W

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